Carony®

Transfer from wheelchair to car seat without lifting.







| SWIVEL BASE COMPATIBILITY | | | | | | | | |
|---------------------------|-----------------|----------------------|-------------------|--------------|------------|--|--|--|
| TURNOUT | TURNY MANUAL | TURNY LOW VEHICLE | TURNY HD/ORBIT | TURNY EVO | 6-WAY BASE | | | |
| × | ~ | × | ✓ 3 | ✓ 3 | × | | | |

3. requires adapter plate with rails

The new Carony, like the old, is a wheelchair that lets the user stay seated in the wheelchair when moving into the car. By docking onto a swivel base the entire wheelchair seat, with the seated user, can slide into the car, where it becomes the car seat.

Our latest version has been improved in many different ways in terms of comfort for the user, ergonomics and functions for the carer as well as a simplified installation process.

The seat itself is highly adaptable with individually adjustable armrests, torso support and thigh support. A built in pneumatic bladder allows the front of the cushion to be raised to further improve seating comfort.

It comes in a discreet grey colour and has a sleek design matching both the interior and exterior of today's cars.

Suitable for people of all sizes, including those of short stature and children.

For a full review of the Carony's features, visit Dealerweb.

Quick facts

- Integrated push handles with brakes
- Many adjustment points for user comfort
- 16" high quality wheels with reflective stripes
- Improved stability
- One model only
- Uses same seat for both wheelchair and vehicle

Specifications

| FRAME MATERIAL | metal |
|----------------|-----------|
| FRAME COLOUR | grey |
| SEAT COLOUR | black |

To guarantee your safety Autoadapt products are designed and tested according to current directives and standards.




Wheel base accessories







Requires a swivel seat for docking to a car.

Recommended accessories: Anti-tip device and Calf/heel support

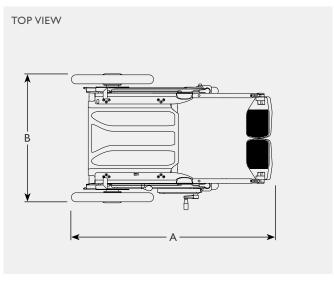

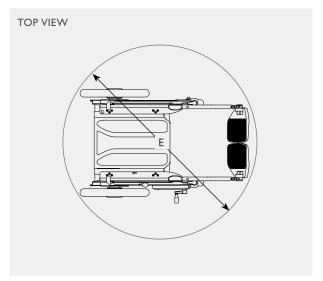

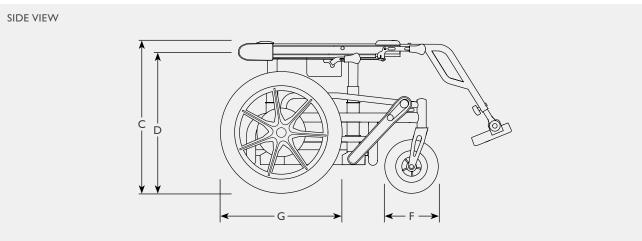

| 104069 | Anti-tip device | 104083 | Thigh supp | ort, right | | 100425 | 2-point hip | belt | | | | |
|--------|-------------------------------------|--------|------------|---------------------|------------------------|----------------------------|----------------------------|---------------------|---------|--------------------------|-----------------|------------------------------|
| 104070 | Puncture free wheels, front | 104084 | Thigh supp | Thigh support, left | | 104056 | 4-point shoulder harness | | | | | |
| 104071 | Puncture free wheels, rear | 104085 | Torso supp | oort, right | : | 104059 | Calf / heel : | Calf / heel support | | | | |
| 104081 | Elevating leg support, right | 104086 | Torso supp | oort, left | | 104057 | Universal el | astic strap | | | | |
| 104082 | Elevating leg support, left | 100665 | Head supp | ort (neck | roll) | 104058 | Elastic exter | nsion strap | | | | |
| Carony | , | | Braye en | Length A | Ard ^{er} B | r _{reight} r C | oo ^{king ke} D | gir Turning E | kront.* | need diameter Rear wh | Assimul Assimul | n leet weight Unit weight |
| 103927 | Carony 16" | | both* | 71 | 64 | 43-57 | 37-51 | 103 | 7" | 16" | 120 | 15.5 |
| 104072 | Adapter plate with rails: Turny Evo |) | | | | | | | | | | |
| 104073 | Adapter plate with rails: Turny HD | /Orbit | | | | | | | | | | |

Seating accessories

Carony® accessories

Getting seated



Anti-tip device

With the intention to make it as easy as possible to use and to adapt the wheelchair to any situation, the anti-tip device is both activated and deactivated using the feet. Be it going down a couple of steps or traversing a higher obstacle.

104069 Anti-tip device



Puncture free wheels

Both front and rear inner tubes are available in a puncture free version. All of them with the same premium quality and look.

104070 Puncture free wheels, front104071 Puncture free wheels, rear



Elevating leg support

Improved leg support for those with long legs or rigid knees.

I04081 Elevating leg support, rightI04082 Elevating leg support, left



2-point hip belt

Provides good support and prevents you from sliding forward in the seat.

100425 2-point hip belt



Thigh and torso support

Fitted anywhere along the slit along the backrest or under the seat, the premium adjustable torso support and thigh support provides a rigid brace wherever they're needed.

I04083 Thigh support, right
I04084 Thigh support, left
I04085 Torso support, right
I04086 Torso support, left



4-point shoulder harness

Distributing pressure across a 25% greater surface area than a traditional H-style harness, the 4-point shoulder harness has a body-contouring shape for optimal support and a comfortable fit. Curved pads follow the contours of the shoulders and ribcage, with a low-placed sternum strap for easy adjustability.

104056 4-point shoulder harness



Universal elastic strap

The 148 cm strap provides temporary or all-day support for travel, work, therapy, sports and fittings. Use on the chest, abdomen, thighs or lower legs. Adds 72 cm to the original Universal elastic strap.

I04057 Universal elastic strap (148 cm) I04058 Elastic extension strap (72 cm)



Calf / heel support

With great pressure distribution and durability, they support the lower legs and control flexion without risk of bruising. Breathable fabric provide air circulation and comfort. Attaches around the leg rests using Velcro. Easy to move.

104059 Calf / heel support

Carony® Classic

A wheelchair that transforms into a car seat, while the user stays seated.







| SWIVEL BASE COMPATIBILITY | | | | | | | | |
|---------------------------|-----------------|----------------------|-------------------|-----------------------|------------|--|--|--|
| TURNOUT | TURNY MANUAL | TURNY LOW VEHICLE | TURNY HD/ORBIT | TURNY EVO | 6-WAY BASE | | | |
| ✓ ² | × | × | ✓ ² | ~ ² | × | | | |

2. requires Carony rails

The Carony Classic is a wheelchair that lets the user stay seated while transferring to the car seat. By docking the wheelchair to either a swivel seat or seat lift, the wheelchair seat can be slid off its base and onto the docking plate where the seat locks safely in place and becomes the car seat. In short, a user can move from wheelchair to car seat and vice versa completely without lifting.

Available with 12" or 24" wheels and suitable for people of all sizes including those of short stature and children.

Quick facts

- Uses same seat for both wheelchair and vehicle
- Highly customisable
- 12" or 24" rear wheels
- Approved to be used as a wheelchair during transport (requires transport fittings)
- Compatible with Tilda tilt function

Specifications

| FRAME MATERIAL | meta |
|----------------|------|
| FRAME COLOUR | red |

To guarantee your safety Autoadapt products are designed and tested according to current directives and standards.


NOTE! See the Carony Classic configurator for available options.







Carony® Classic Dimensions and models

Requires a swivel seat for docking to a car. Recommended accessories: Tilda

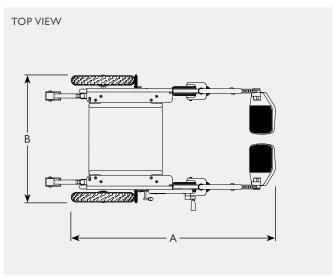

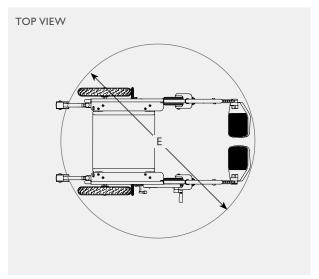

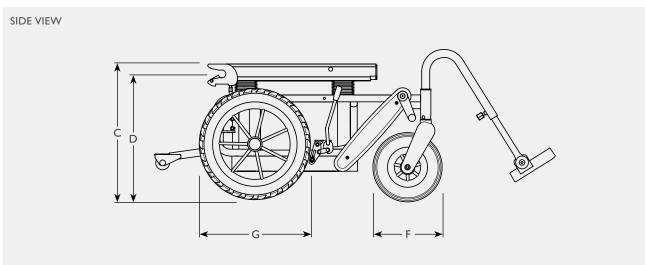

| | Brake lev | erside Length | "Aidely | , keidit | Dodung height | Turning | radius contain | Rear who | ed diameter Matimus | Jules weight |
|-----------------------------------|-----------|------------------|---------|------------|---------------|---------|----------------|----------|------------------------|--------------|
| Carony Classic, chassis | φ. | A | В | C | D | E | F | G | 4. | V |
| Carony Classic 12", chassis | right | 67 | 61 | 38-48 | 35-45 | 84 | 8" | 12" | 100 | 15 |
| Carony Classic 12", chassis | left | 67 | 61 | 38-48 | 35-45 | 84 | 8" | 12" | 100 | 15 |
| Carony Classic 24", chassis | right | 79* | 67 | 38-48 | 35-45 | 84 | 7" | 24" | 100 | 19 |
| Carony Classic 24", chassis | left | 79* | 67 | 38-48 | 35-45 | 84 | 7" | 24" | 100 | 19 |
| Carony Classic 12" fixed, chassis | right | 67 | 61 | 38, 43, 48 | 35, 40, 45 | 84 | 8" | 12" | 100 | 14 |
| Carony Classic 12" fixed, chassis | left | 67 | 61 | 38, 43, 48 | 35, 40, 45 | 84 | 8" | 12" | 100 | 14 |

Carony® Classic Configurator

Choose a model, seat and add docking rails.







| Caron | Class | ssic Fix | ked 12" |
|-------|-------|----------|---------|
| Caron | , Cia | | CG IZ |

INCLUDES

- Seat (choose which)
- Chassis 12" fixed height
- 12" rear wheels
- 8" front wheels
- footrest
- carriage
- push handle
- armrests

Carony Classic 12"

INCLUDES

- Seat (choose which)
- Chassis 12" adjustable height
- 12" rear wheels
- 8" front wheels
- footrest
- carriage
- push handle
- armrests

Carony Classic 24"

INCLUDES

- Seat (choose which)
- Chassis 24" adjustable height
- 24" rear wheels
- 7" puncture free front wheels
- footrest
- carriage
- push handle
- armrests
- anti-tip protection

| SEAT | RIGHT | LEFT | RIGHT | LEFT | RIGHT | LEFT |
|---------------------|--------|--------|--------|--------|--------|--------|
| BEV seat | 101232 | 101233 | 100399 | 100402 | 100872 | 100889 |
| BEV Kids seat | 101245 | 101246 | 101243 | 101244 | 101423 | 101422 |
| Compact seat | 102715 | 102716 | 102711 | 102712 | 102713 | 102714 |
| Recaro Specialist L | 102717 | 102718 | 100683 | 100684 | 100905 | 100906 |
| GS seat | 105027 | 105028 | 105023 | 105024 | 105025 | 105026 |











BEV seat

BEV Kids seat

Compact seat

Recaro Specialist L

GS seat

• Select docking rails for the swivel seat to use it with the Carony Classic

| | Carony Classic Fixed 12" | Carony Classic 12" | Carony Classic 24" |
|----------------|--------------------------------------|--------------------------------------|--------------------------------------|
| SWIVEL BASE | CARONY DOCKING RAILS | CARONY DOCKING RAILS | CARONY DOCKING RAILS |
| Turnout | 475 mm (102720) | 475 mm (102720) | 550 mm (102721) / 600 mm (102722) |
| Turny HD/Orbit | 475 mm (102720) | 475 mm (102720) | 475 mm (102720) |
| Turny Evo | 550 mm incl. adapter plate* (102811) | 550 mm incl. adapter plate* (102811) | 550 mm incl. adapter plate* (102811) |

^{*} On the Turny Evo the Carony rails are mounted on an adapter plate.

Carony® Classic accessories

Getting seated



2-point hip belt

Provides good support and prevents you from sliding forward in the seat.

100425 2-point hip belt



4-point shoulder harness

Posture support in all directions. Curved pads with a low-placed sternum strap for easy adjustability.

104056 4-point harness



Universal elastic strap

Versatile support strap that can be used on the chest, abdomen, thighs or lower legs.

104057 Universal elastic strap 104058 Elastic extension strap



Neck roll

Extra support for the neck and head.

100665 Neck roll



Elevating leg support

Improved leg support for those with long legs or rigid knees.

100424 Elevating leg support, right 100426 Elevating leg support, left



Calf / heel support

Add support to the calfs or place behind the heels to prevent the feet from slipping off the foot rests.

104059 Calf/heel support



Foot support, small

A shorter foot rest for children or people of short stature.

100754 Foot support, small



Heel belt

Adds supports to the heels and prevents the feet from slipping off the foot rests.

100630 Heel belt



Puncture free wheels

The original Carony Classic wheels but filled with gel instead of air.

100438 12" rear, 2pcs (CC 12") **100697** 8" front, Ipcs (CC 12"/24") **101278** 7" front, Ipcs (CC 24")



Anti-tip device

Prevents the Carony from tipping over backwards.

102719 Anti-tip protection



Transport fittings

Tie-down securement points for public/private transport.

101166 Transport fittings



Level raising kit

Raise the seat height 50mm using this kit.

100700 Raising kit, 50mm

• See also seat specific accessories for each seat and the Careva and Crossit posture belts.