

# Carony® Go

A powered wheelchair that transforms into a car seat, while you stay seated.



The Carony Go is a powered wheelchair that lets you stay seated while transferring to the car seat. By docking the wheelchair to either a swivel seat or seat lift, the wheelchair seat can be slid off its base and onto the docking plate where the seat locks safely in place and becomes the car seat. This can be very useful if you are a heavy person and your carer's strength is limited or if moving from the wheelchair is painful for you.

Suitable for people of all sizes including those of short stature and children, the Carony Go has many of the features found on standard powered wheelchairs such as adjustable seat height and narrow wheelbase for indoor use.

### Quick facts

- Uses same seat for both wheelchair and vehicle
- Variable speed setting
- Powered height adjustment
- Signal horn
- Optional head lights and direction indicator lights
- Compatible with Tilda tilt function
- Anti-tip as standard

### Specifications

MAX SPEED .....	6.5 km/h
RANGE .....	25 km
TIME TO FULL CHARGE .....	up to 8 h
MAINS VOLTAGE .....	230 V AC, 50 Hz

To guarantee your safety Autoadapt products are designed and tested according to current directives and standards.



SWIVEL BASE COMPATIBILITY					
TURNOUT	TURNY MANUAL	TURNY LOW VEHICLE	TURNY HD/ORBIT	TURNY EVO	6-WAY BASE
✓ <sup>2</sup>	✗	✗	✓ <sup>2</sup>	✓ <sup>2</sup>	✗

2. requires Carony rails





SWIVEL SEAT	CARONY DOCKING RAILS
Turnout	475 mm (102720)
Turny HD/Orbit	475 mm (102720)
Turny Evo	550 mm incl. adapter plate** (102811)



# Carony® Go

## Dimensions and models

Requires a swivel seat for docking to a car.  
Accessories are available separately.

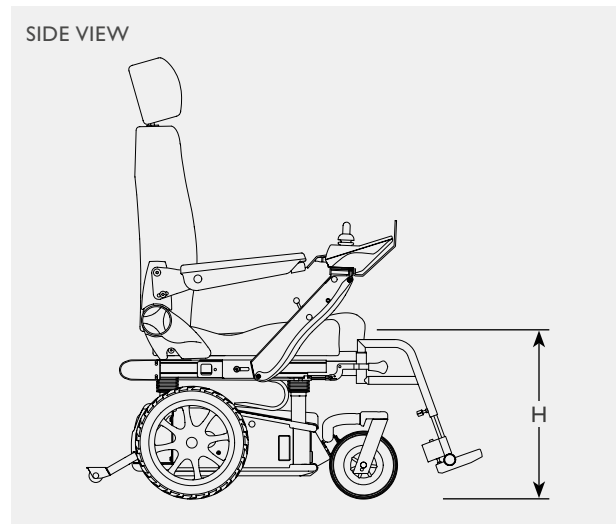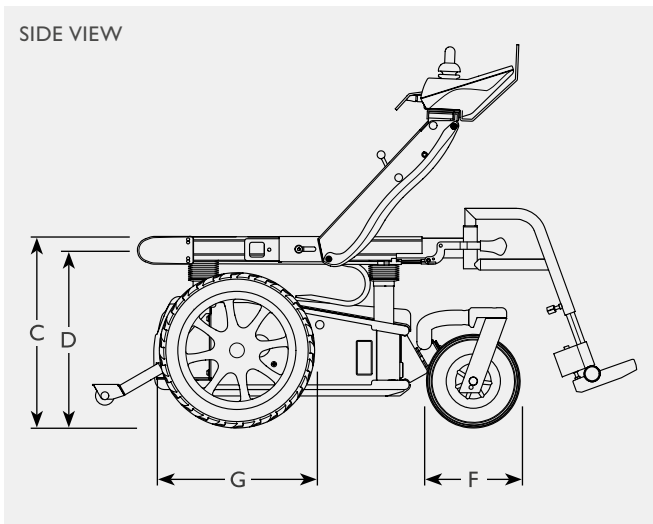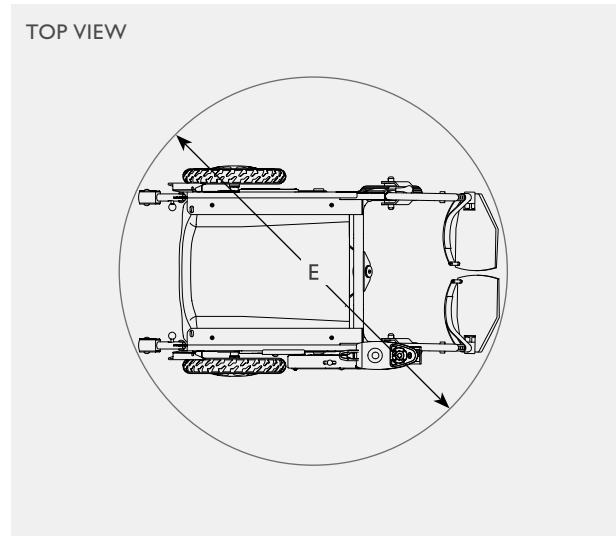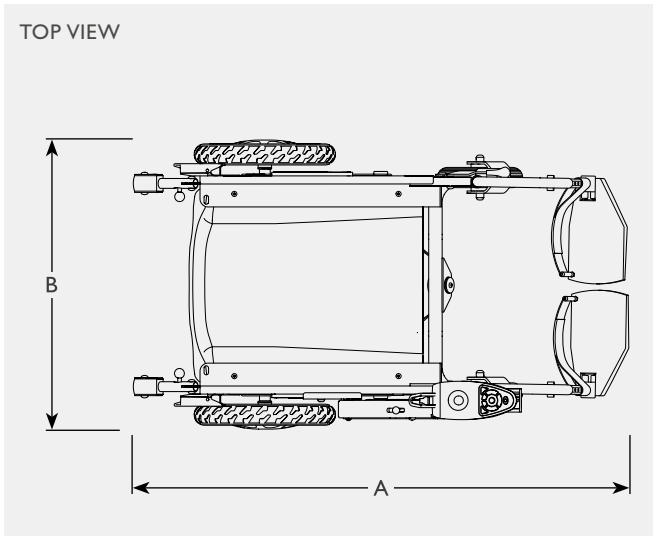

Carony Go		Joystick side	Obstacle height	Length	Width	Height	Docking height	Turning radius	Front-wheel diameter	Rear wheel diameter	Seat height	Maximum user weight	Unit weight
				A	B	C	D	E	F	G	H		
101527	Carony Go	right	8	75	65	38-50	34-46	90	20	31.5	43-55	120	93 / 75
101528	Carony Go	left	8	75	65	38-50	34-46	90	20	31.5	43-55	120	93 / 75
102634	Carony Go, including lights	right	8	75	65	38-50	34-46	90	20	31.5	43-55	120	93 / 75
102635	Carony Go, including lights	left	8	75	65	38-50	34-46	90	20	31.5	43-55	120	93 / 75

All weights and dimensions are given in kilograms and centimetres, and for reference only. Weight capacity refers to the total weight of the seat, accessories and user.

# Carony® Go accessories

## Getting seated

### Carony Go



**Push handle**  
For manual operation.

**101519** Push handle



**Attendance drive**  
A controller for a carer.  
Requires the push handle.

**101542** Attendance drive



**Batteries**  
Replacement batteries (pair).  
12 volt, 45 Ah.

**101520** Carony Go batteries



**Battery charger**  
In: 1x230 volt, 50 Hz, 2 A  
Out: 24 volt, 8 A  
External fuse: 16 A

**101521** Battery charger



**Dynamic programming unit**  
For programming the Carony Go.

**101525** Dynamic programming unit



**Elevating leg support**  
For long legs or rigid knees.

**101575** Elevating leg support, R  
**101576** Elevating leg support, L



**Calf / heel support**  
Add support to the calves or place  
behind the heels to prevent the feet  
from slipping off the foot rests.

**104059** Calf/heel support



**Heel belt**  
Adds supports to the heels and  
prevents the feet from slipping off  
the foot rests.

**100630** Heel belt



**2-point hip belt**  
Provides good support and  
prevents you from sliding forward  
in the seat.

**100425** 2-point hip belt



**4-point shoulder harness**  
Posture support in all directions.  
Curved pads with a low-placed  
sternum strap for easy adjustability.

**104056** 4-point harness



**Universal elastic strap**  
Versatile support strap that can be  
used on the chest, abdomen, thighs  
or lower legs.

**104057** Universal elastic strap  
**104058** Elastic extension strap